SENSITEC NF-3101

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*Congratulations on your purchase of the new Thermometer **Sensitec NF-3101!***

GENERAL DESCRIPTION

The device is designed to measure human body temperature. It is ideal for measuring body temperature of children, even during sleep. You can also measure the temperature of different surfaces, water in the baby bath, baby food and other objects.

We recommend you to know your normal individual body temperature when you feel well in order to properly interpret the data that differ from the average.

Warning: In case if the body temperature is above 38 °C, we recommend you to contact your physician immediately.

I. FEATURES

- Precise non-contact measurement
- Backlit LCD display
- Selectable temperature measurement mode – «**BODY**» or «**SURFACE**»
- User selectable °C or °F
- High temperature alarm
- Memorization of the last 32 measurements
- Automatic data saving and power off

2. SAFETY INFORMATION

- The device must be used only for the purposes described in this user manual
- Do not drop or knock the device, and do not use it in case of damage
- The device must be used only in an ambient temperature range between 10,0 and 40,0 °C
- Do not expose the thermometer to electric shocks
- Do not expose this thermometer to extreme temperature conditions (>50,0 °C or <0 °C)
- Do not use the device in relative humidity higher than 85%
- Do not immerse the device into water or other liquids
- Do not use the device near sources of a high electromagnetic radiation, such as cordless or cellular phones
- Keep the device away from water and heat sources, including direct sunlight
- Do not attempt to charge the batteries that cannot be charged
- Keep the device in the room pending 15-20 minutes before the first use

Before taking of the temperature make sure to remove hair and sweat from the forehead.

Select «**BODY**» mode to measure the body temperature and «**SURFACE**» mode to measure the surface temperature.

Depending on the type and thickness of skin there can be difference in the readings of body temperature.

Should a problem occur with your device please contact your retailer. Do not attempt to repair the device yourself.

3. RECOMMENDATIONS FOR USE

Thermometer Sensitec NF-3101 is designed for non-contact body temperature measurement, namely for forehead and temporal region of adults and children (mode «**BODY**»).

Thermometer Sensitec NF-3101 can also be used to measure the temperature of a baby-bottle or bath, or room temperature (by using the «**SURFACE**» mode).

Thermometer Sensitec NF-3101 is not intended for rectal, oral and axillary method of temperature measurement.

The body temperature of a healthy person can be influenced by the following factors:

- Individual peculiarities of metabolism
- Age (body temperature of infants and small children is higher and it decreases with age. Children are likely to have significant temperature fluctuations faster and more often)
- Clothes that a person is wearing
- Ambient temperature
- Time of day (In the morning the body temperature is lower, and by the end of the day it rises)
- Preceding physical activity
- Measurement method
- Phase of the menstrual cycle


NORMAL BODY TEMPERATURE ACCORDING TO AGE

| Age | Normal temperature, °C |
|---------------|------------------------|
| 0 - 2 years | 36,4 - 38,0 |
| 3 - 10 years | 36,1 - 37,8 |
| 11 - 65 years | 35,9 - 37,6 |
| >65 years | 35,8 - 37,5 |

To achieve the correct and stable results during the first use of the thermometer, we advise to test it and, if necessary, adjust the device with the help of a mercury thermometer (part 8 of this manual).


4. DESCRIPTION

1. IR-Sensor
2. LCD display
3. «MODE SELECTION»
4. «DOWN» button
5. «UP» button
6. Button «MODE»
7. Trigger (power «ON» and measurement button)
8. Battery compartment

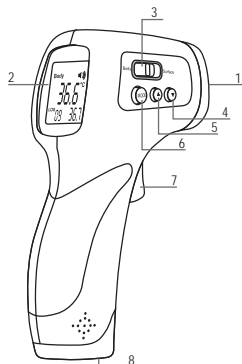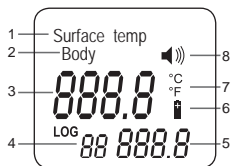5. DISPLAY

1. «SURFACE» mode - surface measurement mode
2. «BODY» mode - body measurement mode
3. Temperature readout
4. The order number of measurement
5. Saved data readout
6. Battery symbol
7. Temperature scale (Celsius of Fahrenheit)
8. Alarm symbol

6. DESCRIPTION OF SYMBOLS ON THE THERMOMETER

| | |
|---|---|
|  | The device is in accordance with Directive EC 93/42/ECC (Medical equipment) |
|  | The device complies with the standard FCC part 15 subpart B: 2007/ Radio Frequency Devices 1C Regulation ICEC-003: 2004 Interference-causing Standard Digital Equipment |
|  | Supply voltage: 3V direct current |
|  | Type of equipment: B (with higher rate of protection) |
|  | With the aim of environment protection, please, utilize the batteries according to the set of rules |
| Body Surface | Identification of measuring method: « BODY », « SURFACE » |
|  | Attention! Read the accompanying documentation |

7. TECHNICAL INFORMATION

| Parameter | Value |
|--|-------------------------------|
| Accuracy of measurement, °C | 0,1 |
| Conditions for use: temperature, °C | 10,0 - 40,0 |
| Storage temperature, °C | 0 - 50,0 |
| Conditions for use: humidity, % | 85 |
| Batteries | DC 3V (2xAA batteries) |
| Size (L x B x H), mm | 149 x 43 x 77 |
| Weight, g | Brutto 400 / Netto 172 |
| Measuring range in « BODY » mode, °C | 32,0 - 42,5 |
| Measuring range in « SURFACE » mode, °C | 0 - 60,0 |
| Measuring distance, cm | 5 - 15 |
| Automatic power off | After 7 seconds of standstill |

PRECISION OF NON-CONTACT MEASUREMENT

| | | |
|----------------------|---------|--|
| From 32,0 to 35,9 °C | ±0,3 °C | According to standard ASTM E1965-1998 (2003) |
| From 36,0 to 39,0 °C | ±0,2 °C | |
| From 39,0 to 42,0 °C | ±0,3 °C | |

8. ADVICE FOR THE FIRST USE ON PRECISION TESTING

Step 1. Take the temperature of a person using a mercury thermometer. You will get 36,5 °C for instance.

Step 2. Compare this result with the normal body temperature, depending on the measurement method as shown in the table below:

| Measurement method | Normal Temperature, °C |
|--------------------|------------------------|
| Rectal* | 36,3 – 37,8 |
| Oral | 36,0 – 37,5 |
| Axillary | 34,7 – 37,3 |
| Ear | 36,0 – 37,8 |

**It should be kept in mind that the rectal method of measurement most closely reflects the central body temperature.*

Step 3. Take the temperature of the same person using the non-contact thermometer keeping the 5 to 15 cm distance between the thermometer and the forehead. Remove any obstacles from the patient's forehead (hair, perspiration, make-up etc). If you get 36,5°C, the non-contact thermometer is properly set and ready for use. If you get a different temperature, 36,4°C (your difference is 1,1°C) for example, you should tune the non-contact thermometer.

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
Step 4. To check, take the temperature again using the non-contact thermometer.

9. PROCEDURE FOR MEASURING BODY TEMPERATURE

9.1. Install the batteries.

9.2. For the first use or when replacing the batteries the device must be stored in the room at least for 15-20 minutes.

- 9.3. If the device has not been used for a long time, we recommend you to measure the room temperature first.
- 9.4. Make sure there is no hair, perspiration, cosmetic or cap covering the forehead.. Point the device at the patient's forehead (see the picture) at a distance of 5 to 15 cm Press the trigger – the temperature will be displayed less than in 1 second.
- 9.5. If there are obstacles for measuring in the forehead area or the ambient temperature is significantly different, the temperature can be measured behind the earlobe.

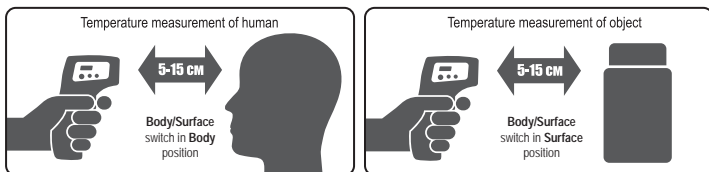10. DIFFERENT MEASURING FUNCTIONS

- 10.1. Choosing The Temperature Unit - «F1» Function.
Press «**MODE**» button for 2 seconds, the screen displays «F1». Select «**DOWN**» for degrees Celsius, «**UP**» for degrees Fahrenheit.
- 10.2. Alarm Setup - «F2».
Press «**MODE**» button for 2 seconds, the screen displays «F1». Press «**MODE**» button again to get «F2». Select «**UP**» to increase the threshold by 0,1°C, and use «**DOWN**» button to reduce it by 0,1°C. By default the signal appears at temperature 38,0 °C.
- 10.3. Adjusting the difference of measurements - «F3» (this function is available only for «**BODY**» mode).
Press «**MODE**» button for 2 seconds, the screen displays: «F1». Press «**MODE**» button two times to get «F3». Select «**UP**» to increase the difference by 0,1°C, «**DOWN**» to reduce it by 0,1°C. Be sure to carry out a verification and adjustment in case of seasonal or environmental changes.
- 10.4. Alarm ON/OFF - «F4» («F3» – for «**SURFACE**» mode).
Press «**MODE**» button for 2 seconds, the screen displays: «F1». Press «**MODE**» button again to get «F4» («F3» – for «**SURFACE**» mode). Select «**UP**» to switch on the buzzer, select «**DOWN**» to switch it off.
- 10.5. Exit the Setting Mode.

Press «**MODE**» button until the screen turns off.

10.6. Data Storage.

All the temperature measurements are saved automatically in the data base, which will display at the right corner of LCD. Press «**UP**» or «**DOWN**» button to display the last temperature measurements. If the thermometer is switched off press «**UP**» and «**DOWN**» buttons together and the last temperature measurement will be displayed. To delete all data chose the memory unit 0 and press the «**MODE**» button.

10.7. Replacement of batteries.

When the LCD displays the battery sign, the battery level is low. Open the lid and change the batteries, pay attention to the correct positioning. Never use rechargeable batteries, use only batteries for single use. If the thermometer is not used for a long time, remove the batteries.

10.8. Longevity of Use.

The non-contact infrared thermometer was designed for professional use, its longevity is guaranteed for 10000 measurements.

II. MAINTENANCE AND CLEANING

- The Infrared Sensor is the most fragile part and must be protected carefully.
- The device can be cleaned with liquid containing 70% alcohol, do not clean the device with corrosive detergent.
- Store the device in a dry environment, and keep it away from water, dust and direct sunlight.

12. TROUBLESHOOTING

If you happen to have one of the following problems, please, refer to this user manual to help resolve the problem. If the problem persists, please, contact our customer service.

- The LCD displays the temperature lower than 32,0 °C. You've probably incorrect switched on the «**SURFACE**» mode, and the thermometer shows the temperature.
- The LCD displays «**HI**». It means that the body temperature is higher than 42,5 °C and the surface temperature is higher than 60,0 °C.
- The LCD displays «**LO**». It means that the body temperature is lower than 32,0 °C and the surface temperature is lower than 0 °C.

13. WARRANTY

- The warranty period for the device is 12 months since the date of purchase. In case if date of purchase is unavailable warranty period starts since the date of manufacturing (see the packaging)
- The warranty is valid only if there is a cashier's receipt or completely filled-in warranty service coupon with the seal of trade enterprise or service centre
- If a manufacturing defect is found during the warranty period the defective unit will be repaired, and if it cannot be repaired – the device will be replaced free of charge
- Manufacturer guarantees compliance with technical characteristics if operating conditions are met
- Warranty service is not performed:
 - when the traces of mechanical action, nicks, cracks, chips, etc. are found on the device (including display and buttons);
 - when the traces of case tampering are found or attempts to repair the device outside the authorized service are disclosed;
 - when the traces of liquid inside the device body are found or if the device has been exposed to aggressive environment;
 - in other cases of violation of the rules of storage, transportation and maintenance of the device set out in the operating instructions.
- The warranty does not apply to the batteries and packaging.
- If during the warranty period, the device has been used for other purposes or the rules of the operating instruction have been violated, repairs will be made at the expense of the consumer.
- If the device was used and can be repaired, it cannot be replaced by a new one.

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Shenzhen Everbest Machinery Industry Co., Ltd., KHP, 19 th